


# Transfer sheet PROFI

*Model 2020*

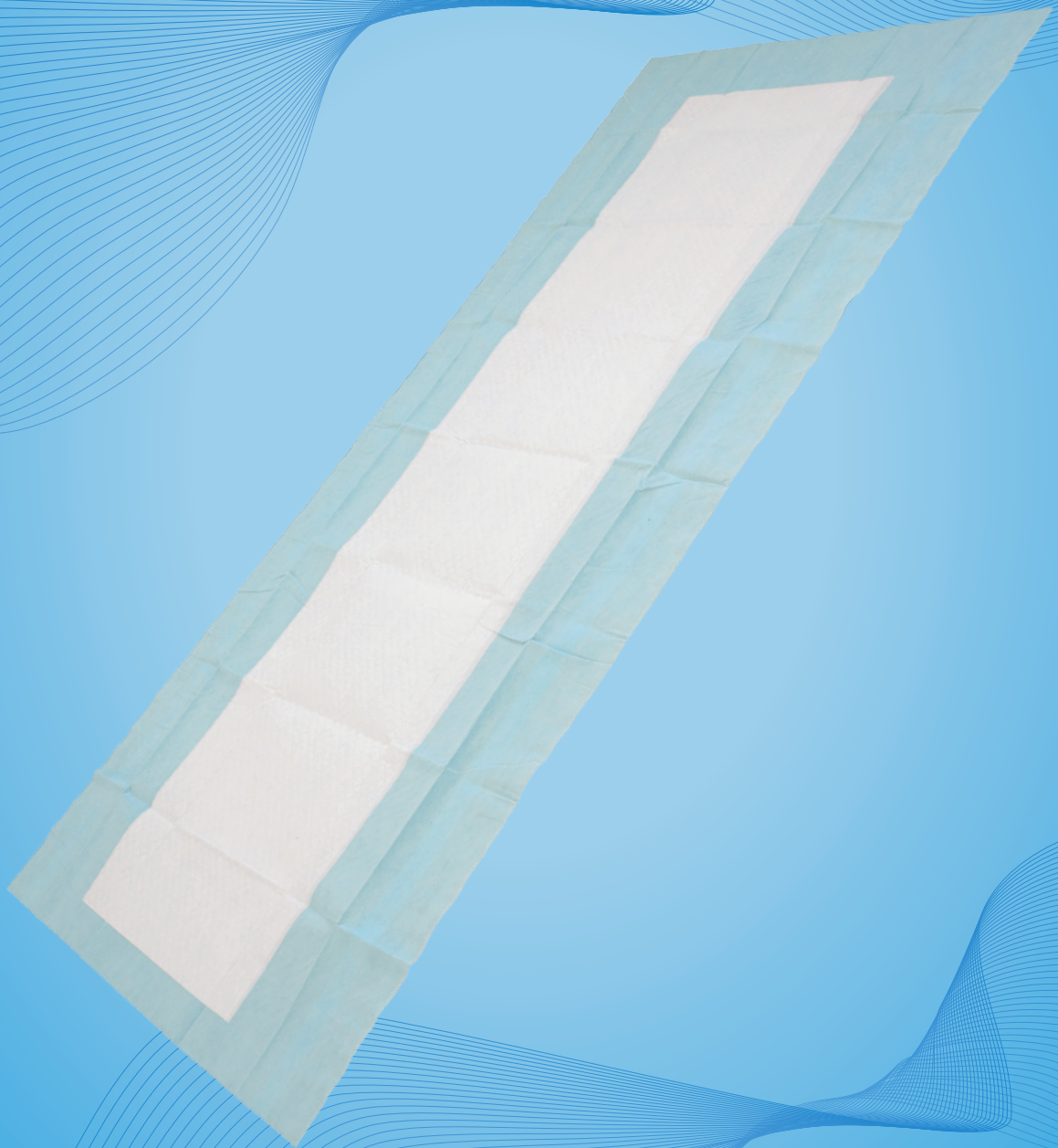

## Extremely resilient and safe

Our innovative HERZmed PROFI transfer sheet is a waterproof single-use sheet for transferring patients for use in rescue services as well as in hospitals and nursing homes.

The sheet is extremely robust, even after absorbing liquid. The integrated, high-performance absorption zone (three layers consisting of PE films, cellulose fluff and fleece) prevents unwanted waterlogging - to the benefit of the patient and the facility.

## Details



Suction area



Backside



Single packaging



## Highlights

### THREE-LAYER SYSTEM:

- 1. PE-film: This material prevents liquids from leaking through. The roughened blue side avoids slipping of the cloth.
- 2. Cellulose fluff: Due to the special rhombic-shaped arrangement, the escaping liquid is quickly distributed in the absorption body.
- 3. Fleece: The upper side consists of a skin-friendly fleece fabric. The fleece layer forms a comfortable patient support with a pleasant feel.

## Extremely resilient

The transfer sheet has a load capacity up to 200 kg.

## High liquid absorption

The transfer sheet has an absorption volume up to 5 liters.

## Cost optimization

Due to the compact packaging design, the transfer sheet takes up little storage space.

## Einzelverpackung

Die Einzelverpackung ermöglicht gute Lagerung, leichte und sichere Handhabung sowie ein hygienisches Arbeiten.

## Fleece upper side

The top side consists of a skin-friendly fleece fabric. The fleece layer forms a comfortable patient support with a pleasant feel.

## Cellulose pulp


The rhombic-shaped cellulose pulp causes an even distribution of the escaping liquid.

## PE-film underside

This material prevents liquids from leaking through and is at the same time responsible for the high load-bearing capacity of the cloth. The roughened blue side prevents the cloth from slipping and thus actively supports the user in his work.

## Three-layer system

The interaction of the three different layers ensures optimal functioning.

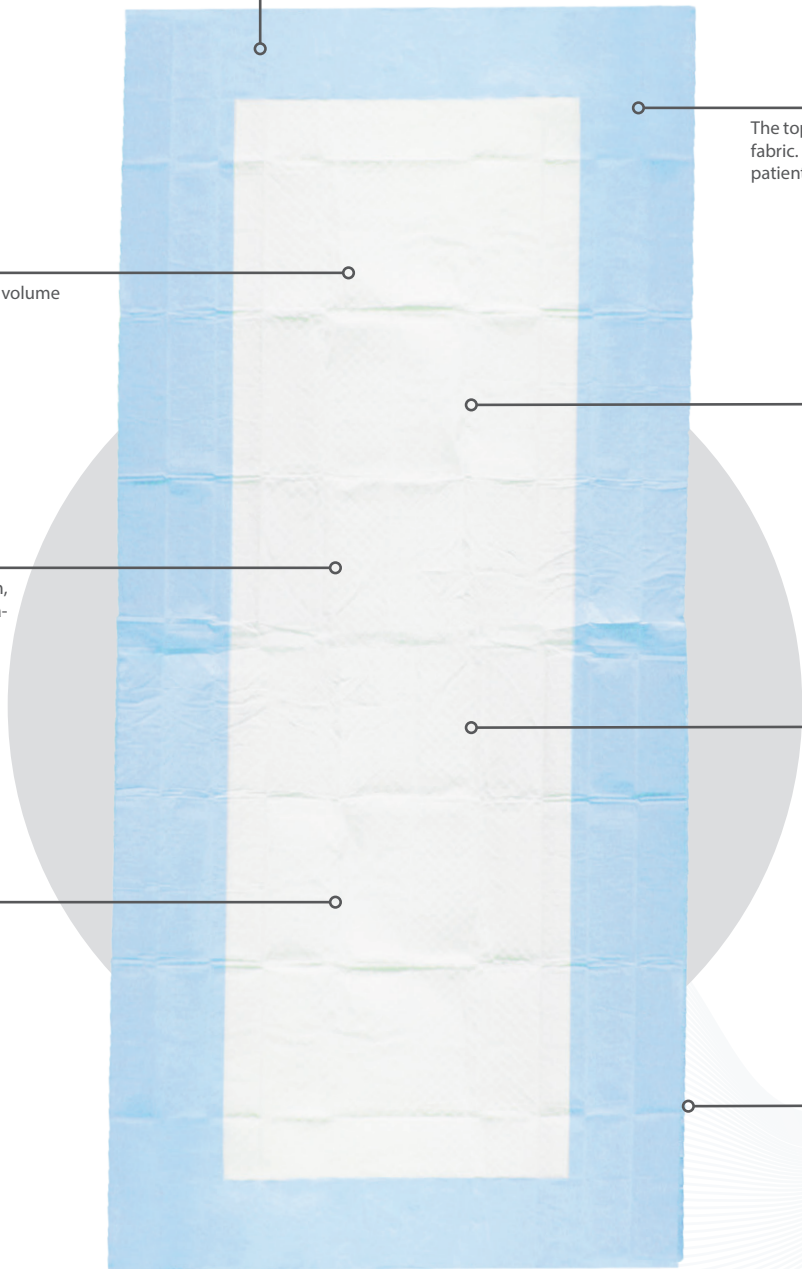

## Hygienic working environment

Due to the generous dimensions, the entire work area is safely covered and effectively protected against damaging influences caused by moisture. The transfer sheet PROFIL can minimize so-called cross-infections for patients and staff.

This is a single-use product!

## Technical specifications

**Size:**

Main fabric: 220 x 110 cm (L x B)  
Suction zone: 200 x 50 cm (L x B)

**Weight:**

Ca. 350 g

**Coloring:**

White/Blue

**Packaging:**

Single packed

**Absorption volume:**

Absorption: Up to 5 liters

**Load capacity:**

Up to 200 kg  
(even after liquid absorption)

**Material:**

Under layer: Polyethylene (PE)  
Middle layer: Cellulose fluff  
Patient layer: Fleece

**Sterility:**

No

## Order information

Art.-No.	Description
V100048	Transfer sheet PROFi, Model 2020

## Order and advice by:

HERZmed GmbH

Herkulan-Schwaiger-Gasse 16  
82487 Oberammergau

**Tel.:** +49 (0) 911 / 14 88 959 - 0

**Fax:** +49 (0) 911 / 14 88 959 - 9

**Mail:** [herzmed@herzmed.de](mailto:herzmed@herzmed.de)

**Web:** [www.herzmed.de](http://www.herzmed.de)

The products shown are only an excerpt of our product range. Exact details and product information are available on request or at [www.herzmed.de](http://www.herzmed.de). All prices marked with an asterisk are exclusive of VAT. All brand names used remain the property of the respective brand owner. We assume no liability for errors, possible misprints and color deviations. Our general terms and conditions apply, which you can also view at [www.herzmed.de](http://www.herzmed.de).

